

OPERATING INSTRUCTIONS

To Start your **NIGHT MAGIC** massage system:

1. Pressing the **ON/OFF** button automatically selects 15 minute timer.
2. Pressing the **TIMER** button switches between the 15 minute and 30 minute timers as indicated by the LED light.
3. Pressing the **AUTO** button (AUTO LED light on) pre-program sets the NIGHT MAGIC massage system to automatically cycle through each of the Massage Therapy Modes for however long the timer has been set.
4. Pressing the **MANUAL** button (AUTO LED light off) allows you to customize your massage session.
5. Pressing the **ZONES** button allows you to select any of 4 zones to target your massage session from head to toes.
6. Pressing the **INTENSITY** button enables you to increase or decrease the intensity of your massage session.
7. Pressing the **MODE** button allows you to select from 9 different pre-programmed Massage Therapy Modes including; Relax, Sleep, Stimulate, Wave and Pulse among the many modes that can be varied with 3 intensity levels and distributed to 4 individual zones as well as all zones so you can customize your massage session to your personal preferences.



NightMagic^{T.M.}

SLEEP MASSAGE THERAPY

by



OWNERS MANUAL OPERATING & INSTALLATION INSTRUCTIONS

Model 2010 Massage System

*LIKE HAVING YOUR OWN
PERSONAL MASSEUSE!*

- Use on Twin, Full, Queen, or King Size Mattress*
- Use on Inner Spring, Foam, or Air Mattress
- Low Voltage Power Supply is UL listed

*2 units recommended with split box spring King or Queen beds.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

DANGER: To reduce risk of electric shock always unplug the unit from the electrical outlet immediately after using and before cleaning.


WARNING: To reduce risk of burns, fire, electric shock, or injury to persons **READ ALL INSTRUCTIONS BEFORE USING.**

- **IMPORTANT:** Patients using a pacemaker should consult a physician prior to use.
- The Night Magic Massager should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Close supervision is necessary when this unit is used by, on, or near children, invalids, or disabled persons.
- Use this unit only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- To disconnect, turn all controls to the off position, then remove plug from outlet.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Not recommended for use with waterbeds.
- For any questions call Customer Support at 1-800-466-6629.

SAVE THESE INSTRUCTIONS

© 2009 by INNOMAX CORPORATION
530 West Elk Place
Denver, Colorado 80216
FAX 303-292-5846
www.innomax.com

INSTALLATION INSTRUCTIONS

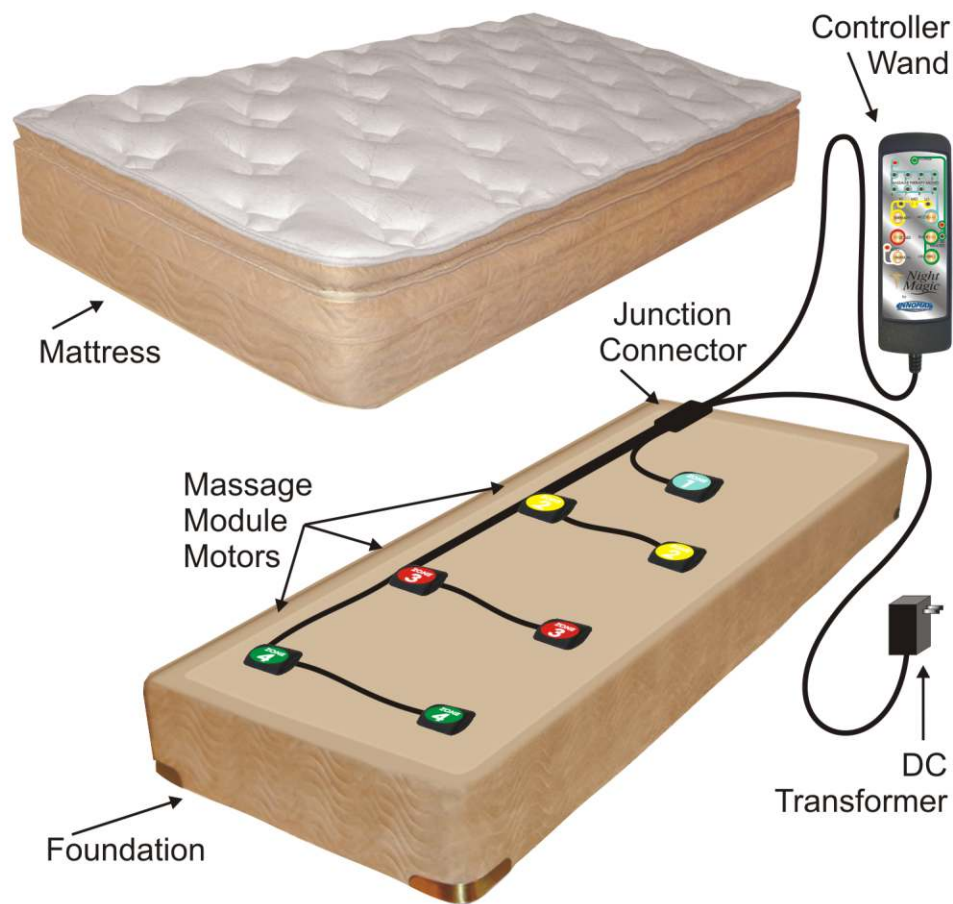

Installation Instructions: For Household Use

1. Do Not Plug Night Magic Massage System into Power Supply until all components have been connected together.
2. Arrange Massage Module Motors on top of the mattress foundation from head to foot with Zone 1 under the neck area, Zone 2 under the shoulders and so forth.
2. Plug Massage Motors into the appropriate Junction Connector.
3. Plug Night Magic Controller Wand into the appropriate Junction Connector.
4. Replace Mattress on top of Foundation.
5. Plug DC Transformer into Power Supply and enjoy your Night Magic Massage System.